This new generation pipette controller fits any hand in a very comfortable way. Shape is ideal for working with glass or plastic pipettes within laminar flow or microbiological safety cabinet. A maximal working efficiency, even with large pipetting volumes, is reached through a powerful, yet guiet motor pump and the availability of several setting options. One-year warranty.

The Profiller™ electro 447

Product advantages

- Powerful pump fills 25 mL in 3 seconds
- Simple, fast operation
- Comfortable finger triggers
- Large display always visible
- Real time battery charge indicator
- NiMH battery rechargeable during use
- Reliable construction for shock and chemical resistance
- Interchangeable protection membrane filter
- Compact table and wall stands



Pipettes up to 100 mL and over







Choice of dispensing modes ①

Instant activation of dispensing with blow-out (TC) or gravity mode (TD).


Instant speed setting ①

A switch enables selection of fast or slow aspiration and dispensing speeds, clearly indicated on display.

Ordering information - Instrument

Supplied with power supply and cable, two spare hydrophobic filters (0.45 µm and 0.2 µm), bench stand, wall hanger and instruction manual.

Profiller™ <i>electro</i> 447	Packaging	Cat. No.
With Europe style plug	1 / pk	447.100E
With UK style plug	1 / pk	447.100G
With USA/Japan style plug	1 / nk	447 100U

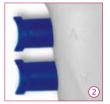

Finger triggers ②

Aspiration and dispensing controlled through two ergonomically positioned finger triggers requiring minimum activation effort. Triggers enable gradual speed regulation.



Battery charge ③

Instrument plug enables charging battery either in storage position or during use. Charge level visible on display.

Convenient stands 4

Instrument supplied with stable bench stand 4 including 3 suckers, and adhesive wall hanger (5).

Ordering information - Accessories and spare parts

Description	Packaging	Cat. No.
Spare bench stand	1 / pk	320.947
Spare wall hanger	1 / pk	320.447
Spare 1.2V NiMH battery	3 / pk	900.917
Power supply, Europe style plug	1 / pk	900.947E
Power supply, UK style plug	1 / pk	900.947G
Power supply, USA-Japan style plug	1 / pk	900.947U

Description	Packaging	Cat. No.
Membrane filter, 0.45 µm	5 / pk	322.447
Membrane filter, 0.2 µm	5 / pk	322.407
Silicone pipette holder	1 / pk	1.447.01
Nose cone	1 / nk	1 447 02

